



platinum  
STAIRLIFTS

CURVE

custom built comfort



[www.platinumstairlifts.com](http://www.platinumstairlifts.com)

# CURVE

tailor made for  
a smooth ride

## curve standard

- Ultra thin, twin tube steel rail
- Quiet and smooth ride quality
- Ergonomically shaped joystick control
- Tight internal bend on a 200mm radius
- Fits tight to the staircase
- Revolutionary ERGO chair
- Digital diagnostic display in arm
- Remote controls
- Retractable reel seat belt
- Suitable for users up to 120kg (19st)

## curve optional features

- Powered swivel and footrest options
- Lever linked footplate option
- Foot cover option
- RF remote control
- Bariatric Curve option available up to 160kg (25st)
- ERGO Space seat for narrow staircases, saving 50mm with a standard ERGO seat pad or 75mm with an ERGO Plus seat pad.



Ergonomic Joystick




Chair Swivel



V4 Adjustable Seat Plate



Digital Diagnostics

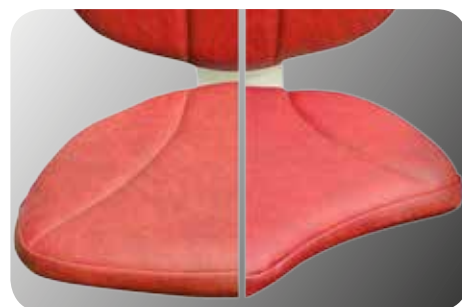

ERGO Seat Pad



# custom built comfort

The Platinum Curve is hand built in the UK, to the exact measurements of your staircase, with the carriage being fully adjustable to your size and body shape.

Twists, turns and corners are no problem for your Curve, as the rail is tailor made to your individual staircase. With the Platinum Curve, you get a smooth, safe, comfortable ride and the freedom of your house once more.

In addition, the Curve's compact design and floor-mounted rail means it takes up as little space as possible and does not affect your home's decor.

For curved staircases, there is no better option than the Platinum Curve.



ERGO chair with  
Space seat and  
ERGO Plus pad.  
See back page  
for more details.



# CURVE

reduce a whole staircase  
down to five simple steps



1

Lower the arms and the seat pad, put the footplate down and your Curve stairlift is ready to use. The digital diagnostic screen will display `Platinum` meaning you are ready to go.





2

Your Curve chair will be fully adjusted to your size and shape at installation, so all you need to do is strap yourself in, sit back, relax and enjoy the ride!




3

To operate your Platinum Curve stairlift, simply use the ergonomic joystick in the arm and lightly press in the direction you wish to travel.



4

At the top of the stairs, use the side-mounted swivel lever (or optional powered swivel), to turn the chair and safely step out onto the landing.

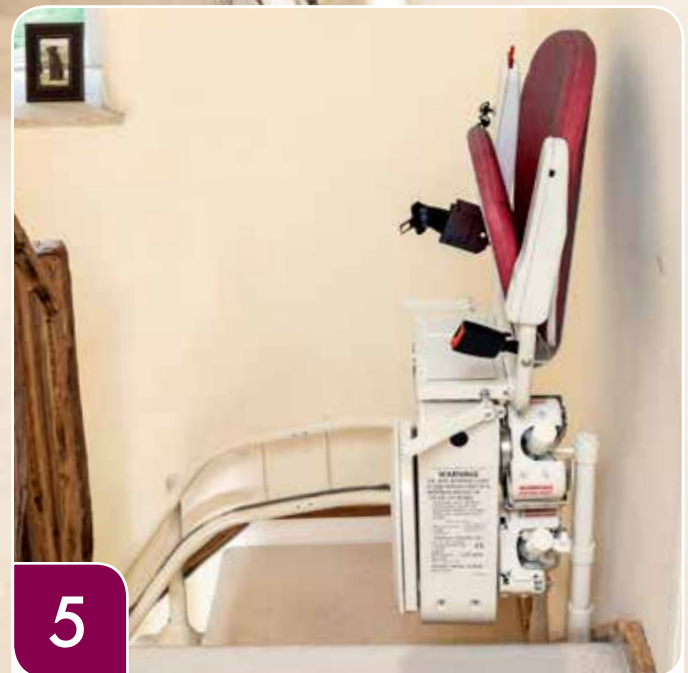

5

When not in use, your Curve stairlift can be folded and parked on the stairs, leaving clear access for other family members.

# ask for ergo

The Platinum ERGO seat is one of the most advanced stairlift seats available. The ERGO is fully adjustable to the size and shape of the user and comes with a wealth of built-in features. The ERGO seat comes as standard on the Curve. Ask your distributor for more details.



**Ergonomic Joystick**  
Total control at your fingertips.



**Security Lock**  
Prevent unauthorised use of the lift.



**ERGO Plus Pad**  
Maximum comfort for shorter legs.



**V4 Adjustable Seat Plate**  
Allows for horizontal adjustment of the seat back. Combined with the multi-position seat pad, the ERGO offers multiple configuration options.



**Digital Arm Display**  
Easy to read, real-time diagnostics.



**Retractable Safety Belt**  
Additional safety when lift is in operation.



**ERGO Standard Pad**  
Designed for maximum comfort.



**Lever Linked Footrest**  
Opens and closes footrest without the need to bend down. Powered option also available.



**Safety Edge**  
Senses obstructions and brings the chair to a stop.



# one seat hundreds of configurations

The ERGO is the only stairlift seat that can be fully adjusted to suit the user at point of installation. When installed with the Curve stairlift, the ERGO offers unparalleled support and comfort for users up to 120kg (19st).

*The ERGO has been designed with over 800 configurations, with multiple positioning for the arms and back rest, together with variations on swivel, foot rest and upholstery.*

Depending on the height of the user, the amount of space needed for the stairlift to fit safely on the staircase can alter by a matter of inches - a significant consideration on narrow stairs. To accommodate this, the ERGO seat has two different seat pad options.

## ergo standard seat



## ergo plus seat



# ergo space less is more



The compact ERGO Space seat sits just 40mm from the wall, providing an extra 50mm of space compared to the standard ERGO seat. When combined with the ERGO Plus seat pad, the Space seat can reduce the area needed to operate by a further 25mm, a total saving of 75mm.

The Space seat's back rest is 80mm shorter than the standard ERGO seat, making the ride more comfortable as users sit further back into the seat and making the whole carriage less visually impactful on decor.

The ERGO Space Seat - a smaller, more comfortable chair, giving you a bigger, more usable space. Less is definitely more.



[Modernstairlifts.com](http://Modernstairlifts.com)  
0863118566

